


## Horton Technologies Adjustable Click Bar

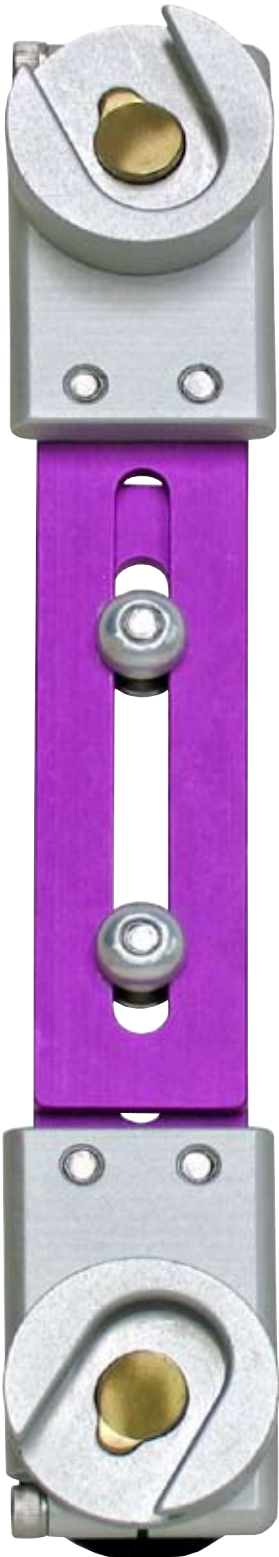

The Horton Click Bar was designed by an orthotist who understood the need for a more user friendly device when treating a child with the Ponseti method.

The Click Bar allows feet to be rotated internally or externally and is easy to use. The shoes are fitted first, then clicked into place on the bar and adjusted to the required angle. This offers excellent foot positioning within the shoe and reduces the chances of the child kicking them off during fitting. Donning and doffing is much easier too; with just a click you can detach the bar from the shoes, change the child then click the shoes back on to the bar.

Available in two sizes and with the option of a 10 degree kicker bar.



Adjustable Click Bar is recommended for use with

**A-Line**  
by Blatchford



6"-8" Adjustable Bar: FR408

8"-12" Adjustable Bar: FR412

6"-8" Adjustable Bar with 10° Kicker: FR408-10

8"-12" Adjustable Bar with 10° Kicker: FR412-10