


ProMax Rheumatoid Insole

The ProMax Rheumatoid Insole offers superior shock absorption, comfort and cushioning. Designed as a pressure redistribution insole and developed from high spec. materials, ProMax helps relieve pain, decrease pressure and reduce joint damage.

Interface Features

- Dual layer construction
- Lunatec motion SL top layer
- Closed cell structure, durable and hygienic
- Absorb shearing forces
- Excellent bedding propertise in the horizontal load plane

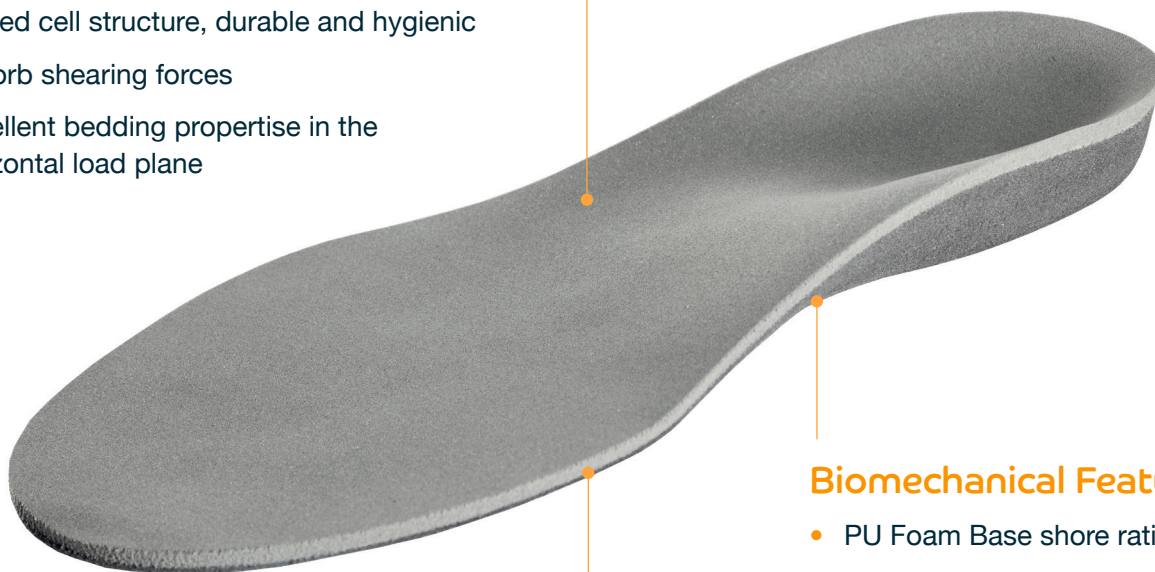

Increased Patient Comfort

- PU base layer assists shock absorption
- Specified PU foam reduces and transfers unwanted pressure from prominent bone structures
- Excellent for Rheumatism or sensitive diabetic feet
- Reduces pain through minimising shear forces for pressure sensitive area

Biomechanical Features

- PU Foam Base shore rating A40-50
- Shock absorbing material not intended for rear foot posting
- Excellent abrasion and tear resistance
- Very high impact resistance while providing comfort and stability for the patient

Full insole additions are also available for selection on RMS including forefoot additions incorporated into the price. Please choose 'ProMax Rheumatoid Insole' in the product section of RMS to specify this insole.