


ProTec III Diabetic Insole

The ProTec III Diabetic Insole is a pressure relieving insole designed to decrease the rate of ulceration, enhance comfort and reduce shear force on the foot. It provides the same benefits of the standard ProTec II Diabetic Insole but features an additional triple layer for enhanced durability.

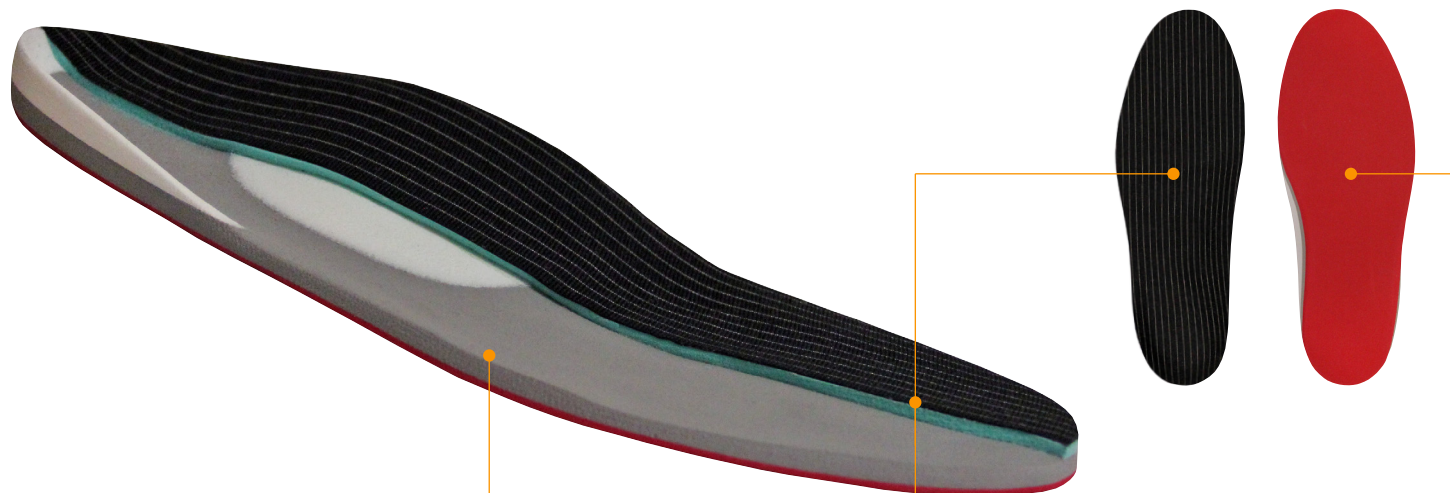

Biomechanical Features

- Trilaminar EVA base, Shore densities are A65-30-50
- High memory Poron 96 mid layer (optional)
- X-STATIC® top layer

Interface Features

- X-STATIC® top cover - antimicrobial fabric clinically proven to reduce HAI-causing contaminants. 99.9% pure metallic silver bonded fibers create an ionic shield, inhibiting the growth of bacteria and fungi. (N.B. not suitable for deep heel cups, for use on shallow trims only)
- Anti-fungal properties
- Thermo regulator

Increased Patient Comfort

- Excellent impact absorption
- Cushioning microcellular structure
- Open cell - breathable

Full insole additions are also available for selection on RMS.

Please choose 'ProTec III Diabetic Insole' in the product section of RMS to specify this insole.